



RESISTENTI
AGLI IMPATTI



ROBUSTI



LAVABILI



ANTIBATTERICI




RESISTENTI RAGGI UV



AMAGNETICI

BRM 002

NON-MAGNETIC PATIENT TROLLEY

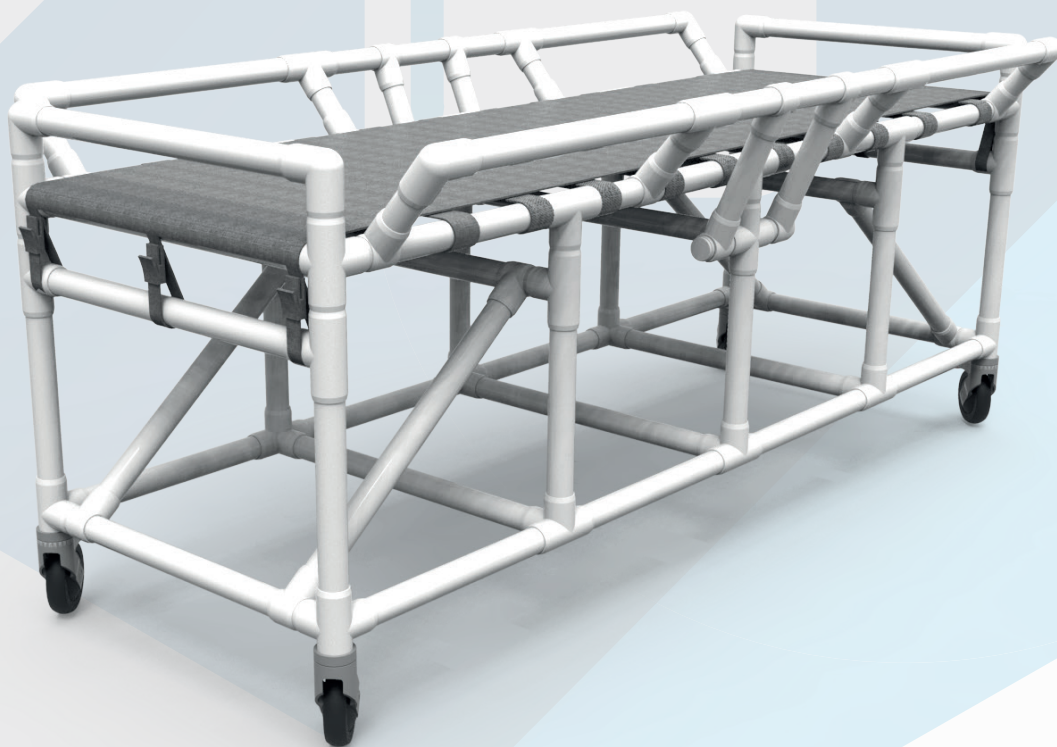

ENCLARA
MEDICAL



The BRM 001 non-magnetic patient trolley is designed for easy transfer of patients with limited functional mobility into restricted access zones of MR environments.

BRM 001 is strictly made of non-metallic, non-magnetic, non-conductive materials. For this reason, it is fully compliant with the requirements set by the EN 62570:2016 harmonized standard practice for MR SAFE marking

FRAME

The frame is entirely made of the best furniture grade PVC. The addition of titanium dioxide in the pvc compound results in a better appearance by giving it a UV-stable, glossy-white finish; it also enhances mechanical resistance and provides specific antibacterial properties.

Components are chemically welded together using a PVC-base glue, so that the frame behaves like a single part, solid, rigid and sturdy.

PATIENT SUPPORT

The patient lies on a comfortable, breathable support, made of hypoallergenic (class I according to the OEKO-TEX standard 100) and fireproof (non-ignition according to the EN 1021-1 and 1021-2 standards) polyester mesh. The support is easily removable for optimal sanitization.

SIDE RAILS

The side-rails can be unlocked and folded down with a simple operation performed while standing in a central position, which allows the assistant to keep control over the patient.

ACCESSORIES



IV Pole



Side rail

